




Iris



Polaris




Isa

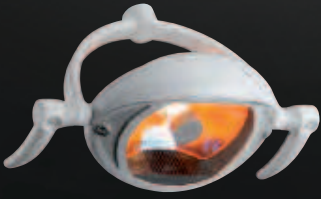

Vision



Nuvola

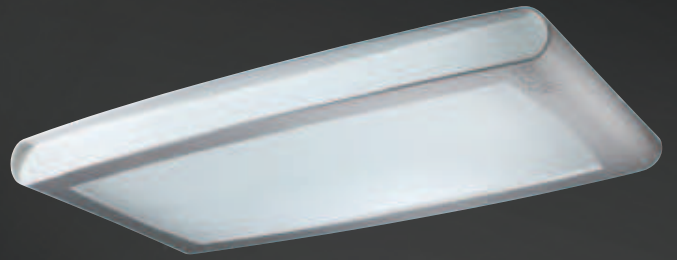

i light it ...

The "I light it" philosophy has revolutionised the way to design and manufacture our products. The goal is to achieve a significant improvement of the dentist's work, paying attention to the quality of light and protection of the environment.

Harmony in lighting is a perfect combination between environmental illumination and the treatment area illumination, our systems could work as one, independently, but in symbiosis to each other. With Gcomm you have now available all the lighting solutions.



Iris

app for innovative communication




Gcomm introduces, first in the dental lighting sector, an innovative user interface with iOS technology, available now for iPad. It represents a revolutionary and unique step towards the traditional dental light concept.

The application "My light" is equipped with an easy interface for the remote regulation of colour temperature, light intensity and Full HD videocamera zooming. There are two more additional options: User customisation and Gcomm recommended values.

Three buttons allow to flip, reverse and freeze the image from the videocamera.

More user profiles can be customised for a maximum flexibility and statistic information.

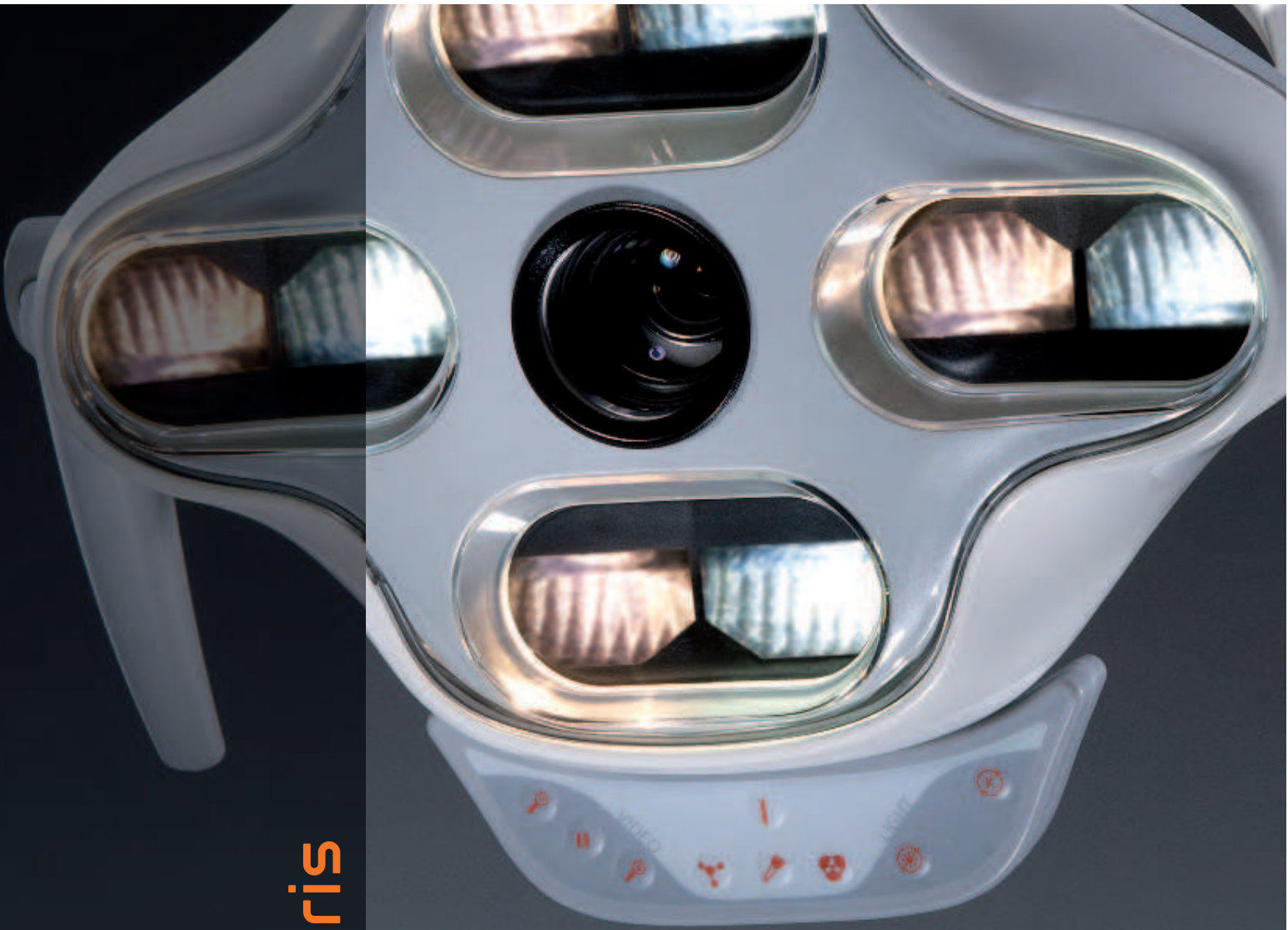

Iris built-in full HD videocamera

IRIS is equipped with a Full HD videocamera (1.920 x 1.080 px) with Autofocus and white balance regulating device. With a 30x optical zoom the videocamera ensures clear and detailed images of the oral cavity for diagnostic and treatment. Videos can be recorded on external devices* for educational purposes (university and conferences) and for the treatment technique improvement.

* third party products




Iris


new performance in led light



Full HD
videocamera system

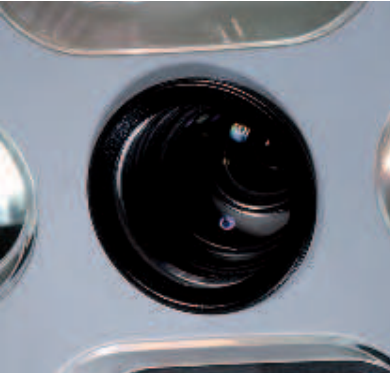

180° keyboard
rotation

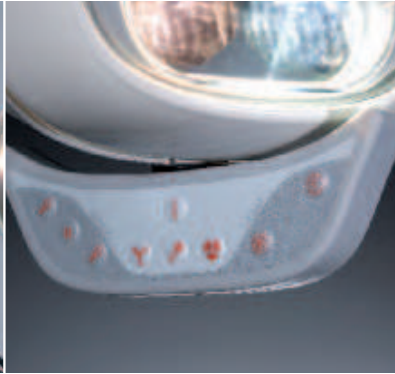

Autoclavable
handles



Patient
mirror



IRIS, the second generation of Gcomm LED light, has multi purpose performances, combining a direct LED illumination technology together with an High Definition (Full HD) video camera system for professional use. Special made lenses together with 8 cold and warm white LEDs ensures a perfect light pattern and uniformity in illuminance distribution. A 180° rotation of the keyboard allows the operator to use it in different positions and makes easier the selection of the 3 pre set programmes:

Surgery treatment

colour temperature setted on 4.500 K for a better contrast on soft tissue

Colour tooth matching

colour temperature setted on 5.500 K for a daylight colour rendering

Polimerysation filter

colour temperature setted on 3.000 K to prevent activation of composites

Adjustable light intensity between 8.000 lux and 35.000 lux. Possibility to adjust as well the colour temperature (from 4.200 K to 6.000 K) to the individual needs of the user. Additional keys on the keyboard allow the use of the Full HD videocamera system by zooming in and out, up to 30 x. The main body of IRIS made in aluminium grants a long durability and robustness. Ergonomic and simple design applies to different surgical requirements quickly. Colour rendering index (CRI) >90 % (from 4.200 K to 5.500 K).

POLARIS R2

| 4 | 3 | 2 | 1 |



led technology to **light** your work



The POLARIS R2 light is realised with an advanced design, characterised by simple and soft lines. Design has been particularly studied in order to realise a dental light not only with high performances but also aesthetically considered: a perfect union of functionality and aesthetics.

Soft lines and smooth surfaces, the possibility of an easy handles' extraction and sterilisation allows to an optimal cleaning and hygiene.

POLARIS R2 is characterised by shells entirely in aluminium, which assures a high robustness. The studied geometry of reflectors, in combination with each LED source, allows to realise a luminous flux which is homogeneous, clean and without shadows.

POLARIS R2 is equipped with two potentiometers which allow to adjust both light intensity (from 8.000 lux to 35.000 lux) and colour temperature (from 4.200 K to 5.800 K).



Potentiometer for Lux and Kelvin



Sensor On/Off



Autoclavable handles



ISA dental halogen light



Dental light ISA has been designed to be ergonomic and modern at the same time. The results of the applied technology are:

- an operating light very easy to handle, stable and unique for its armonious lines and pleasant colours
- a perfect light pattern and cold white light (up to 30.000 lux)
- new high performance resin reflector, easy to clean and resistant to brakes
- a colour rendering index > 90

Dental light ISA is the expression of a product MADE IN ITALY.